






US EPA ARCHIVE DOCUMENT

## BIOCHEMISTRY DEPARTMENT

PAGE 1 of 25	METHOD No. AG-335	SUBJECT DETERMINATION OF METHIDATHION AND SOME OF ITS METABOLITES IN MILK BY GAS CHROMATOGRAPHY
EDITION 1/9/79	Replaces AG-98	
SUBMITTED BY: R. A. Kahrs		APPROVED BY: <i>R.T. Murphy</i>

1.0 SCOPE

This method describes procedures for the analysis of milk for the following compounds:

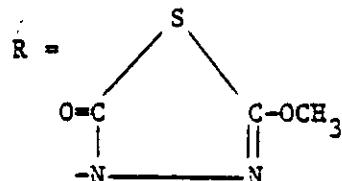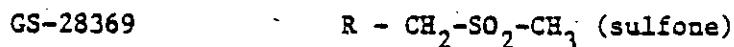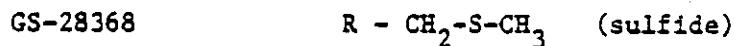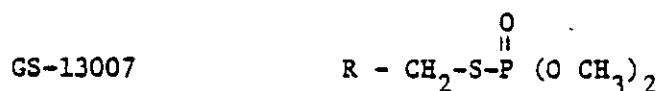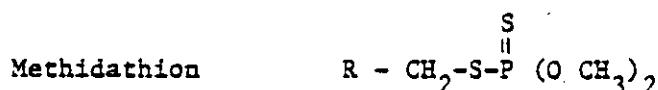2.0 PRINCIPLE OF THE METHOD

Acetone is added to the milk to extract methidathion and its degradation products. The precipitated milk solids are filtered off and these are extracted again with toluene to complete the extraction.

The acetone-water and the toluene extracts are combined in a separatory funnel and shaken. Two phases are formed. The toluene phase (upper phase) contains the methidathion, GS-13007, sulfide, sulfone and most of the sulfoxide. Some of the sulfoxide remains in the aqueous phase. The toluene phase is evaporated to dryness and taken up in n-hexane. This solution is extracted with acetonitrile to recover the compounds sought while the milk fat and other interfering materials remain in the hexane.

## BIOCHEMISTRY DEPARTMENT

PAGE 2 of 25	METHOD No. AG-335	SUBJECT DETERMINATION OF METHIDATHION AND SOME OF ITS METABOLITES IN MILK BY GAS CHROMATOGRAPHY
EDITION 1/9/79	Replaces AG-98	
SUBMITTED BY: R. A. Kahrs		
		APPROVED BY:

The acetonitrile is evaporated to dryness and then the sample is taken up in toluene for gas chromatography to determine parent compound, oxygen analogue, sulfide, sulfone and part of the sulfoxide. The aqueous phase is extracted with methylene chloride to recover the remainder of the sulfoxide. This is evaporated to dryness and taken up in toluene for gas chromatography to determine the remainder of the sulfoxide.

3.0 REAGENTS

Acetone:	Nanograde
Toluene:	Nanograde
Hexane:	Nanograde
Acetonitrile:	Nanograde
Dichloromethane:	Spectroquality
Sodium Sulfate:	Reagent Grade
Methidathion:	Analytical Standard
GS-13007:	Analytical Standard
GS-28368:	Analytical Standard
GS-28369:	Analytical Standard
GS-28370:	Analytical Standard

4.0 EQUIPMENT

Jars:	16 oz flint, wide mouth
Buchner Funnel:	6 cm
Filter paper:	Whatman No. 1
Vacuum adapter:	For erlenmeyer flasks: similar to catalog No. K-20500 adapter, 24/40, Kontes Technical Glassware catalog TG-20, Kontes Glass Co. Vineland, NJ.
Filter vac:	Rubber seal between Buchner funnel and vacuum adapter
Separatory funnels:	125, 250 and 500 ml
Air manifold:	N-EVAP by Organamation or equivalent
Gas Chromatograph:	Microtek MT-220 equipped with a Flame Photometric Detector

## BIOCHEMISTRY DEPARTMENT

GREENSBORO, N.C.

PAGE 3 OF 25	METHOD No. AG-335	SUBJECT DETERMINATION OF METHIDATHION AND SOME OF ITS METABOLITES IN MILK BY GAS CHROMATOGRAPHY
EDITION 1/9/79	Replaces AG-98	
SUBMITTED BY: R. A. Kahrs		APPROVED BY:

5.0 PROCEDURE5.1 PREPARATION OF SAMPLE

Milk samples to be analyzed should be thawed slowly at a temperature between 20-40°C to prevent excessive coagulation and souring. The entire sample should be thoroughly mixed by a vigorous shaking to insure uniformity in sampling.

5.2 EXTRACTION

A 50 gram sample (50 ml) of milk is transferred to a 16 oz wide mouth jar equipped with a plastic cap. A polyethylene liner is placed under the cap to avoid dissolving extraneous material from the cap lid, and to avoid loss of solvent. Two hundred ml of acetone are added to the jar and the mixture is shaken for 1/2 hour. The acetone-aqueous extract is filtered through two pieces of filter paper in a Buchner funnel under suction. The acetone-aqueous extract is then transferred to a 500 ml separatory funnel. The resulting filter pad (not paper) is reshaken with 10 ml of distilled water and 50 ml of acetone for 1/2 hour.

The acetone-water extract is again filtered and combined in the separatory funnel. The resulting filter pad is shaken a final time with 100 ml of toluene for 1/2 hour. The toluene extract is filtered similarly to the acetone extract. The pad is washed with a 75 ml portion of toluene. The toluene and acetone-aqueous extracts are combined and vigorously shaken for three minutes, and the phases allowed to separate.

## BIOCHEMISTRY DEPARTMENT

GREENSBORO, N.C.

PAGE 4 of 25	METHOD No. AG-335	SUBJECT DETERMINATION OF METHIDATHION AND SOME OF ITS METABOLITES IN MILK BY GAS CHROMATOGRAPHY
EDITION 1/9/79	Replaces AG-98	
SUBMITTED BY: R. A. Kahrs		APPROVED BY:

5.3 PARTITION CLEAN-UP5.3.1 Toluene-Acetone Phase:

The toluene-acetone phase is evaporated to an oily residue which is then dissolved in 50 ml of hexane, and transferred to a 125 ml separatory funnel. A 25 ml portion of acetonitrile is used to rinse the flask and is added to the separatory funnel, and the entire mixture is shaken for one minute. The acetonitrile partitioning is repeated twice more with fresh 10 ml portions of acetonitrile. The acetonitrile fractions are combined in a 125 ml separatory funnel and washed with a fresh 50 ml portion of hexane by shaking for one minute. The hexane wash is discarded.

5.3.2 Aqueous Phase:

The aqueous phase is transferred to a 250 ml separatory funnel and partitioned with dichloromethane. Three successive two minute extractions using 50 ml portions of dichloromethane are dried through a 1 inch pad of sodium sulfate. The pad is then rinsed with a 25 ml portion of dichloromethane. This partitioning procedure cannot be too severe or an emulsion will form. The use of an end over end movement in rotating the funnel will minimize the chance of an emulsion.

5.4 DETERMINATION OF METHIDATHION, GS-13007, GS-28368, GS-28369,  
GS-283705.4.1 Designation of Sample Fractions:

In this part of the procedure, there are now two fractions of sample extract to analyze. The acetonitrile fraction will contain methidathion, GS-13007, GS-28368, GS-28369, and 50% of the GS-28370 present in the sample, and therefore, a multiple analysis will be performed. The dichloromethane fraction contains only the remaining 50% of the GS-28370 present in the sample, requiring only a single analysis.

## BIOCHEMISTRY DEPARTMENT

GREENSBORO, N.C.

PAGE 5 of 25	METHOD No. AG-335	SUBJECT DETERMINATION OF METHIDATHION AND SOME OF ITS METABOLITES IN MILK BY GAS CHROMATOGRAPHY
EDITION 1/9/79	Replaces AG-98	
SUBMITTED BY: R. A. Kahrs		APPROVED BY:

5.4.2 Determination Procedures for Both Fractions:

The acetonitrile or dichloromethane is evaporated to 2 ml using a flash evaporator. The residue is transferred quantitatively with dichloromethane to a 2 dram vial, and taken to dryness. The residue is dissolved in an appropriate amount of Nanograde toluene. The final determination is done by means of a gas chromatograph equipped with a Flame Photometric Detector, with a sulfur specific filter (394  $\mu\text{m}$ ). The conditions used for gas chromatography are given in Table I.

The gas chromatographic method is standardized by injecting known amounts of the various compounds. Stock solution of these are made by dissolving 100 mg of compound in 100 ml of acetone. Dilutions are made from this using toluene to obtain solutions ranging from 1 to 3 ng per  $\mu\text{l}$  in concentration. A constant volume of 5  $\mu\text{l}$  equivalent to 5 to 15 ng of compound are used for standardization. Typical chromatograms are shown in Figure 1 and 1-A. The Flame Photometric Detection system, when employed in the sulfur specific mode, has been shown to be volume dependent. In order to obtain uniformity in the response of the detector, only equal volumes of standard and sample are injected (i.e., every injection in the sulfur specific mode is done with 5  $\mu\text{l}$  injections).

A standard graph is prepared from the chromatograms of standards. Peak areas are calculated by multiplying the peak height by the base at half-height. These are plotted against the weight of compound injected. Peak areas of unknown samples are compared directly with such a graph to obtain the amount of compound present in the sample.

## BIOCHEMISTRY DEPARTMENT

GREENSBORO, N.C.	PAGE	METHOD No.	SUBJECT
	6 of 25	AG-335	DETERMINATION OF METHIDATHION AND SOME OF ITS METABOLITES IN MILK BY GAS CHROMATOGRAPHY
	EDITION	Replaces	
	1/9/79	AG-98	
SUBMITTED BY:		R. A. Kahrs	
		APPROVED BY:	

6.0 RECOVERY STUDIES

The recovery of methidathion and metabolites added before extraction to raw cow's milk are shown in Table II, and typical chromatograms shown in Figure 2, 3, 4, 5, 6. This table shows the sample code, the interval in days before and after the start of dosing, the individual compounds that were added, the amount of compound added (in terms of ppm), grams injected, the ppm of each compound detected in the acetonitrile and dichloromethane, the combined ppm found and the percent recovered.

The recoveries from cow's milk are shown at the level of 0.01 ppm of added compound. The pre-treatment samples represent an extraction by blending. No decrease in the percent recovery is obtained when using a shaking extraction. From this data, a lower limit of detectability for raw cow's milk was set at 0.005 ppm for three metabolites and at 0.01 ppm for methidathion and GS-13007 because of their longer elution time. The recoveries of individual compounds are (a) methidathion 88% to 120%, (b) GS-28368 60% to 70%, (c) GS-28369 135% to 160%, (d) GS-28370 130% to 190% of theoretical (refer to section 8.0 Discussion of Method).

7.0 ANALYSIS OF TREATED SAMPLES

For the complete description of sample coding, see the Dawson Report<sup>1</sup>. The TM refers to milk from the cow which was dosed, the CM to milk from the control cow. Minus (-) refers to pre-treatment, and plus (+) to post-treatment. A, refers to AM samples and P, to PM samples, the cows being milked at 6 am and 6 pm. There was one exception, the 18th day of post-treatment (+18). On this day, the treated cow was milked every 6 hours. These milk samples were designated as follows:

TM + 18 F (12 N)	12:00 PM
TM + 18 F ( 6 P)	6:00 PM
TM + 18 F (12 M)	12:00 AM
TM + 18 F ( 6 A)	6:00 AM

The "T" cow was dosed every morning at 8:00 AM.

## BIOCHEMISTRY DEPARTMENT

PAGE 7 of 25	METHOD No. AG-335	SUBJECT DETERMINATION OF METHIDATHION AND SOME OF ITS METABOLITES IN MILK BY GAS CHROMATOGRAPHY
EDITION 1/9/79	Replaces AG-98	
SUBMITTED BY: R. A. Kahrs		
		APPROVED BY:

Raw data obtained in the analysis of treated cow's milk are shown in Table III. This table provides data showing sample code, dosage rate of compound fed to cows in capsule form, the number of days after initial dosing, the number of hours after daily dosage when sample was taken, the weight of sample injected, the individual compounds being analyzed for, the amount (ng) of compound detected in acetonitrile and dichloromethane fraction, the ppm of each compound detected, and the combined ppm. Typical chromatograms are shown in Figure 3.

8.0 DISCUSSION OF METHOD

In setting the lower limit of detectability of 0.005 ppm for the three metabolites, various factors of compound stability were taken into account. When GS-13005 is fed to rats, three metabolites were found by Esser et al.<sup>2</sup> in urine. They were the sulfide (GS-28368), the sulfoxide (GS-28370), and the sulfone (GS-28369). The degradation of methidathion by trans-methylation first produces the sulfide, which being unstable, readily oxidizes to the sulfoxide and subsequently, to the sulfone. Because of the ease of oxidation in milk, recoveries for the sulfoxide and sulfone should be somewhat higher when the milk is fortified with all three standards. As shown in Table II, the recoveries for these compounds are more than 100%. Using the same extraction procedure, but not in the presence of milk, the expected recoveries of 90% to 100% are obtained.

9.0 REFERENCES

1. Murchison, T. E., Dawson Research Corp., "Metabolism Study of  $C^{14}$  Methidathion in a Cow", February 18, 1969.
2. Esser, H. O., Muche, W., and Alt, K. O., Helvetica Chemica Acta, 51, Fasc. 3, 513 (1968).

TABLE I GAS CHROMATOGRAPHIC CONDITIONS

Instrument: Microtek MT 220 equipped with Flame Photometric Detector. Sulfur specific filter (394  $\mu$ ).

Column: 3% SE 30 + 0.3% Epon 1001 on Gas Chrom "Q" (60/80 mesh) packed in Pyrex tubing (2 ft x 1/4 in).

Injector Temperature: 225°C

Column Temperature: 120°C & 170°C

Transfer Temperature: 250°C

Detector Temperature: 160°C

N<sub>2</sub> Carrier: 100 ml/min

O<sub>2</sub> Flow: 30 ml/min

Air Flow: 70 ml/min

H<sub>2</sub> Flow: 150 ml/min

Attenuation: Input: 10<sup>3</sup>, Output: 8

Minimum Detection Limit: Methidathion-----5 ng  
GS-28368-----5 ng  
GS-28369-----5 ng  
GS-28370-----5 ng

Volume Injected: 5  $\mu$ l

Retention Times: Methidathion @ 170°C-----2.5 min.  
~~GS-28368 @ 120°C-----1.1 min.~~  
GS-28369 @ 170°C-----1.1 min.  
GS-28370 @ 170°C-----1.7 min.

Chart Speed: 1/2 inch per min.

Response: Peak areas in inches

NOTE: A 10% DC-200 on Gas Chrom "Q" gas chromatographic column may also be used for these determinations.

AG-335

9 of 25

1/9/79

TABLE II RECOVERIES OF METHIDATHION AND ITS DEGRADATION PRODUCTS ADDED TO MILK

Sample Code	Experiment Day	Compound Determ.	PPM Added	Wt. Inj. (g)		PPM Found		Total PPM	% Recovery
				$\text{CH}_3\text{CN}$	$\text{CH}_2\text{Cl}_2$	$\text{CH}_3\text{CN}$	$\text{CH}_2\text{Cl}_2$		
CM-3 Figs. 2,3,4	- 3	Methidathion	0.00	1.00	2.50	<0.01	<0.01	<0.01	--
			GS-28368	0.00	0.99	--	<0.005	--	<0.005
			GS-28369	0.00	1.00	2.50	<0.005	<0.005	<0.005
			GS-28370	0.00	1.00	2.50	<0.005	<0.005	--
CM-3 plus 0.01 ppm* Figs. 2,3,4	- 3	Methidathion	0.01	1.00	2.50	0.009	<0.01	0.009	90
			GS-28368	0.01	0.78	--	0.007	--	0.007
			GS-28369	0.01	1.00	2.50	0.013	<0.005	0.013
			GS-28370	0.01	1.00	2.50	0.010	0.005	0.015
CM+17A Figs. 4,5,6	+17	Methidathion	0.00	1.20	1.20	<0.01	<0.01	<0.01	--
			GS-28368	0.00	1.15	1.20	<0.005	<0.005	<0.005
			GS-28369	0.00	0.78	1.20	0.006	<0.005	0.006
			GS-28370	0.00	1.20	1.20	<0.005	<0.005	--
CM+17A plus 0.01 ppm* Figs. 4,5,6	+17	Methidathion	0.01	0.60	1.20	0.012	<0.01	0.012	120
			GS-28368	0.01	1.14	1.20	0.006	<0.005	0.006
			GS-28369	0.01	0.60	1.20	0.016	<0.005	0.016**
			GS-28370	0.01	0.60	1.20	0.010	0.006	0.016

\* 0.01 ppm of each of the 4 compounds was added

\*\* When corrected for the check this value becomes 0.010 ppm (100% recovery).

AG-335

10 of 25

1/9/79

TABLE III RAW DATA FROM DETERMINATION OF METHIDATHION AND ITS DEGRADATION PRODUCTS IN MILK

Sample*	Experiment Code	Day	Dosage mg/kg/day	Hours After Dosage	<u>Wt. Inj.</u> (g)	Compound Determined	<u>Ng Found</u>		<u>PPM Found</u> <u>CH<sub>3</sub>CN</u> <u>CH<sub>2</sub>Cl<sub>2</sub></u>	<u>Total PPM</u>
							<u>CH<sub>3</sub>CN</u>	<u>CH<sub>2</sub>Cl<sub>2</sub></u>		
TH-3	- 3D	1	---	1.0	2.5	Methidathion	< 5	< 5	< 0.01	< 0.01
				1.0	-	GS-28368	< 5	--	< 0.005	< 0.005
				1.0	2.5	GS-28369	< 5	< 5	< 0.005	< 0.005
				1.0	2.5	GS-28370	< 5	< 5	< 0.005	< 0.005
TM+17P Figs. 7, 8	+17D	1	12	0.6	1.2	Methidathion	< 5	< 5	< 0.01	< 0.01
				1.1	0.8	GS-28368	< 5	< 5	< 0.005	< 0.005
				0.3	1.2	GS-28369	6.9	5	0.023	0.023
				0.3	0.8	GS-28370	10.3	12.3	0.035	0.035
TM+17A Figs. 7, 8	+17D	1	24	1.2	1.2	Methidathion	< 5	< 5	< 0.01	< 0.01
				1.1	1.2	GS-28368	< 5	< 5	< 0.005	< 0.005
				0.8	1.2	GS-28369	5.7	< 5	0.007	0.007
				1.2	1.2	GS-28370	< 5	5.3	< 0.005	< 0.005(t) < 0.005(t)
TM+18F (12H)	+18D	1	18	0.35	-	Methidathion	< 5	--	< 0.01	< 0.01
				1.2	-	GS-28368	< 5	--	< 0.005	< 0.005
				0.35	-	GS-28369	3.9	--	0.011	0.011
				0.35	-	GS-28370	5.8	--	0.017	0.017
TM+18F (6P)	+18D	1	12	0.30	-	Methidathion	< 5	--	< 0.01	< 0.01
				1.2	-	GS-28368	< 5	--	< 0.005	< 0.005
				0.30	-	GS-28369	4.7	--	0.016	0.016
				0.30	-	GS-28370	4.9	--	0.016	0.016
TM+18F (12N)	+18	1	6	0.35	-	Methidathion	< 5	--	< 0.01	< 0.01
				1.2	-	GS-28368	< 5	--	< 0.005	< 0.005
				0.35	-	GS-28369	9.8	--	0.028	0.028
				0.35	-	GS-28370	9.8	--	0.028	0.028

AC-335



11 of 25

1/9/79

TABLE III RAW DATA FROM DETERMINATION OF METHIDATHION AND ITS DEGRADATION PRODUCTS IN MILK (cont'd)

Sample*	Experiment Code	Day	Dosage mg/kg/day	Hours After Dosage	Wt. Inj. (g)	Compound Determined	NH <sub>3</sub> Found CH <sub>3</sub> CN / CH <sub>2</sub> Cl <sub>2</sub>	PPM Found CH <sub>3</sub> CN / CH <sub>2</sub> Cl <sub>2</sub>	Total PPM
TM+18F (6A)	+18	1	24	0.53	--	Methidathion < 5	< 0.01	--	< 0.01
				1.1	--	GS-28368 < 5	< 0.005	--	< 0.005
				0.53	--	GS-28369 4.9	0.009	--	0.009
				0.57	--	GS-28370 < 5	< 0.005	--	< 0.005
TM+25P Fig. 9	+25	1	12	0.96	1.5	Methidathion < 5	< 0.01	< 0.01	< 0.01
				1.2	0.2	GS-28368 < 5	< 0.005	< 0.005	< 0.005
				0.32	1.5	GS-28369 6.7	0.021	< 0.005	0.021
				0.32	0.9	GS-28370 8.6	9.1	0.027	0.037
TM+28P Figs. 10, 11	+28	2	12	1.2	1.2	Methidathion < 5	< 0.01	--	< 0.01
				0.56	--	GS-28368 4.6	--	0.008	--
				0.19	--	GS-28369 11.4	--	0.061	--
				0.19	--	GS-28370 11.3	--	0.061	--
TM+28A Fig. 10	+28	2	24	1.1	--	Methidathion < 5	< 0.01	--	< 0.01
				0.35	--	GS-28368 4.3	--	0.012	--
				0.22	--	GS-28369 4.7	--	0.022	--
				0.22	--	GS-28370 5.6	--	0.026	--
TM+32P Fig. 12	+32	2	96	0.60	1.25	Methidathion < 5	< 0.01	< 0.01	< 0.005
				1.25	1.25	GS-28368 < 5	< 0.005	< 0.005	< 0.005
				0.60	1.25	GS-28369 2.7	< 5	< 0.005	< 0.005
				0.60	1.25	GS-28370 5.3	5.3	< 0.005	< 0.005

\* See Section 7.0 for Sample Code explanation.

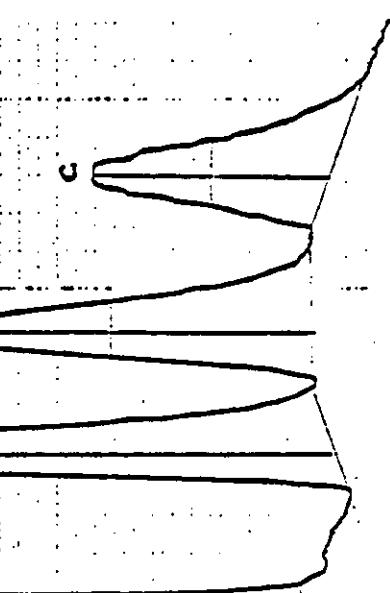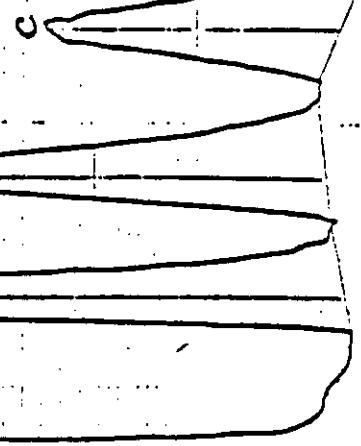

**Fig. 1. Standardization of GLC Response**

- A. GS-28369
  - B. GS-28370
  - C. Methidathion
- Column Temperature: 170°C
- 1 7.5 ng 11 10 ng III 12.5 ng


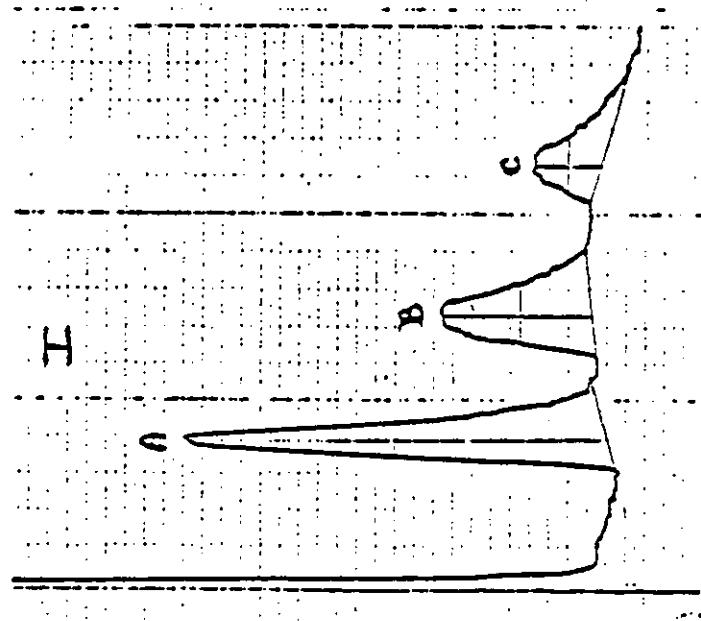


Fig. 1-A  
Standardization of GIG Response to GS-28168 (sulfide).

1 5 ng      II 7.5 ng      Column Temperature: 120°C

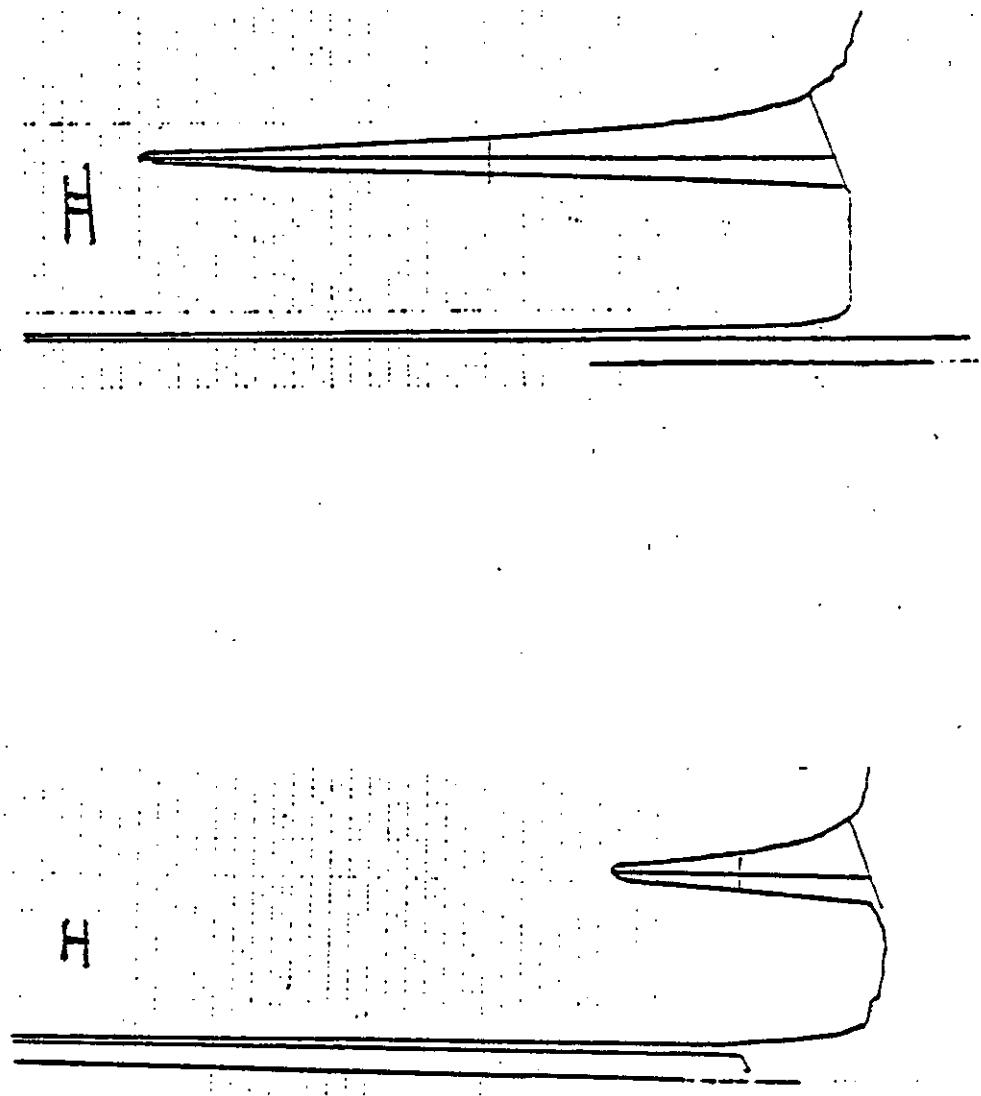

AG-315

13 of 25

1/9/79

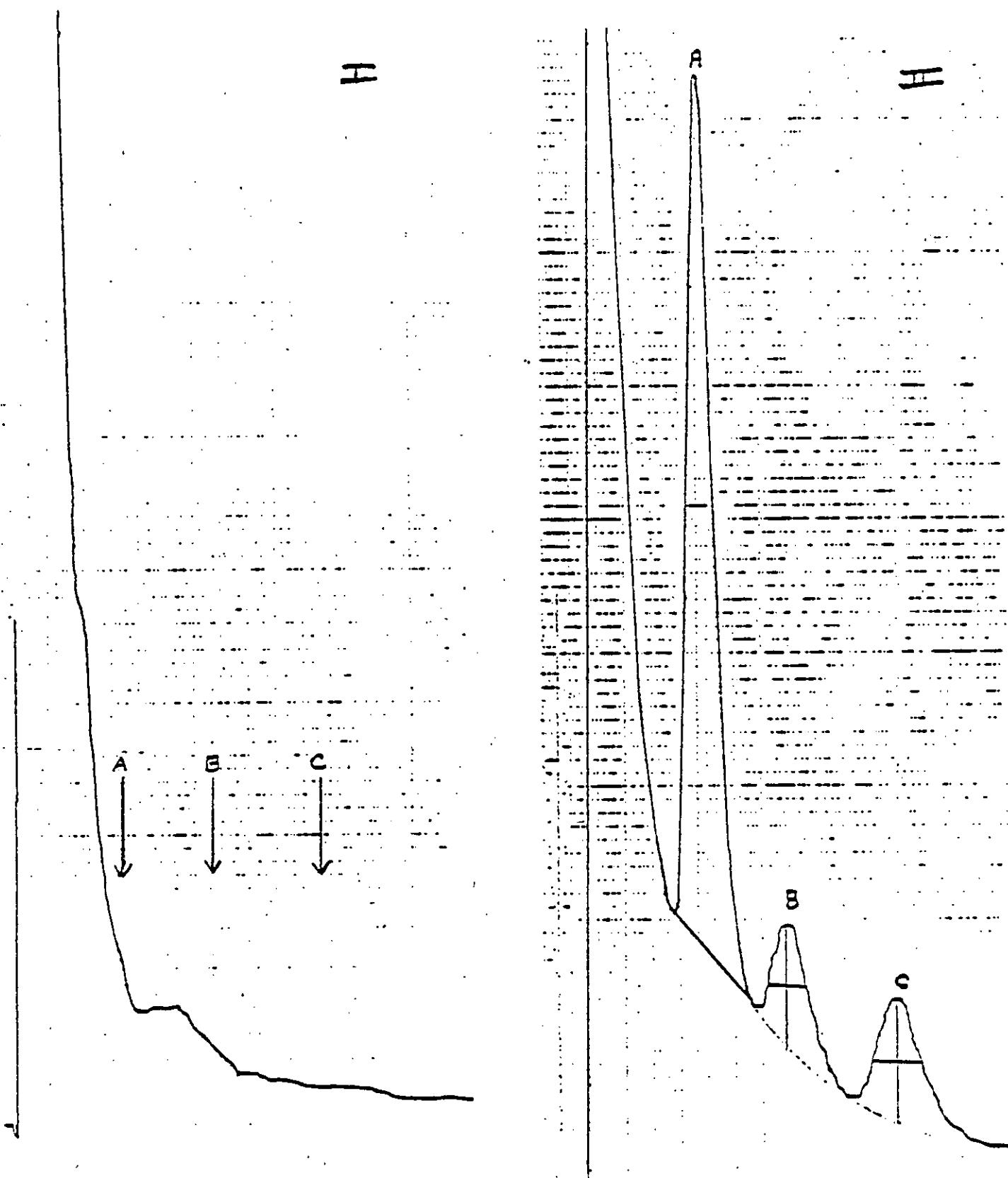


Fig. 2. Gas Chromatograms of Recoveries of Added Compounds (Acetonitrile Fraction).  
A. GS-28369      B. GS-28370      C. Methidathion  
I. Check Milk CM-3  
II. Check Milk CM-3 plus 0.01 ppm of each compound.

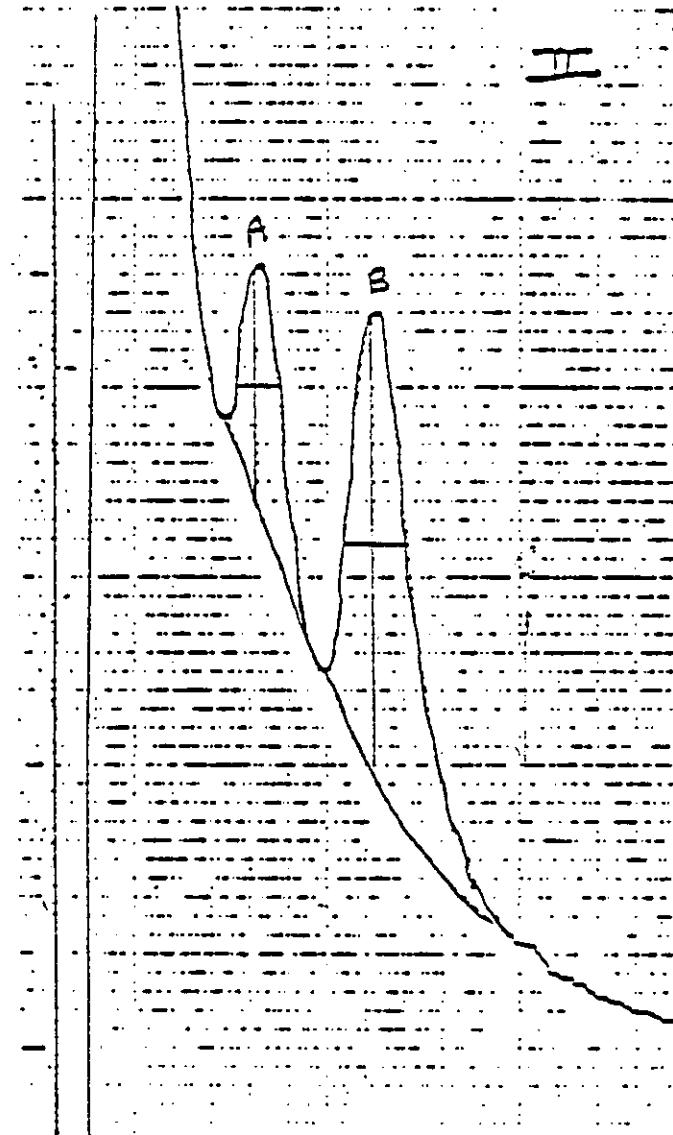
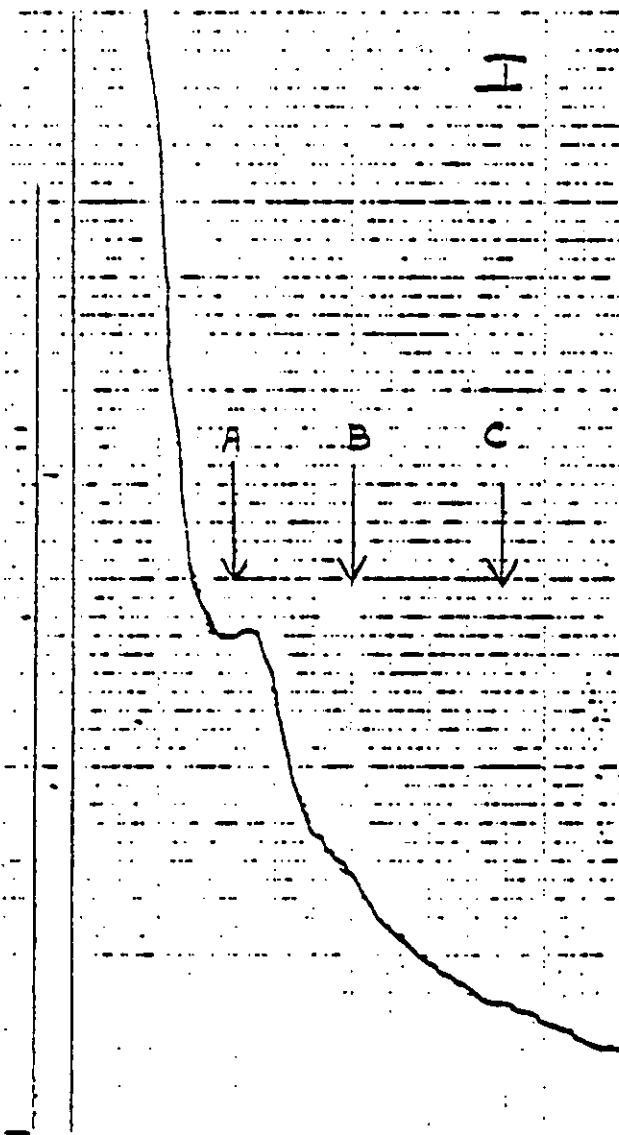


Fig. 3

Gas Chromatograms of Recoveries of Compounds Added to Milk (Methylene chloride fraction).

A. GS-28369

B. GS-28370

C. Methidathion

I. Check milk Q4-3

II. Check milk Q4-3 plus 0.01 ppm of each compound.

AG-335

16 of 25

1/9/79



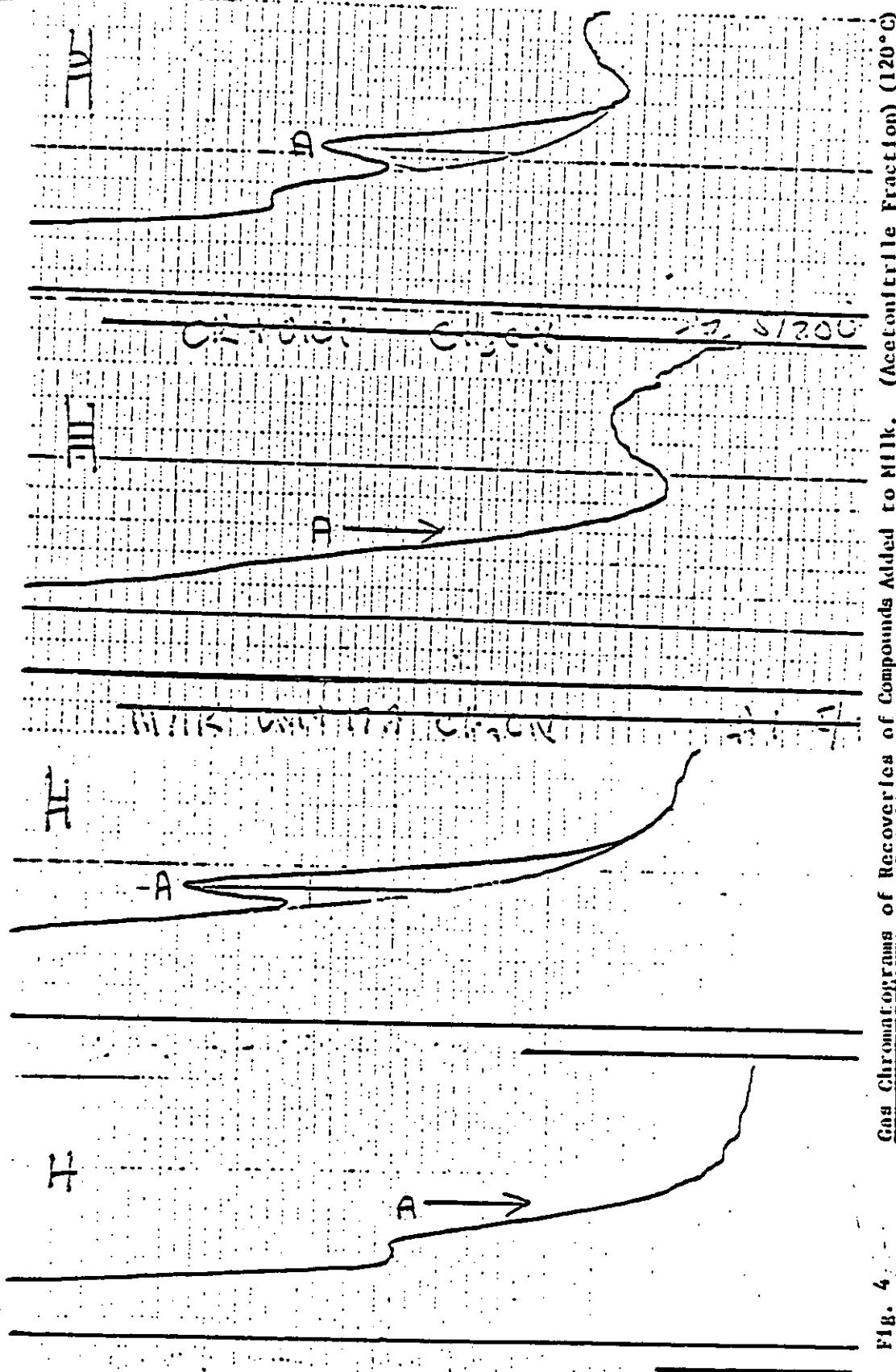


Fig. 4 Gas Chromatograms of Recoveries of Compounds Added to Milk. (Acetonitrile Fraction) (120°C)

- I Check Milk CM-3
- II Check Milk CM-3 plus 0.01 ppm of GS-2836B (D)
- III Check Milk CM-17A
- IV Check milk CM-17 plus 0.01 ppm of GS-2836B (D)

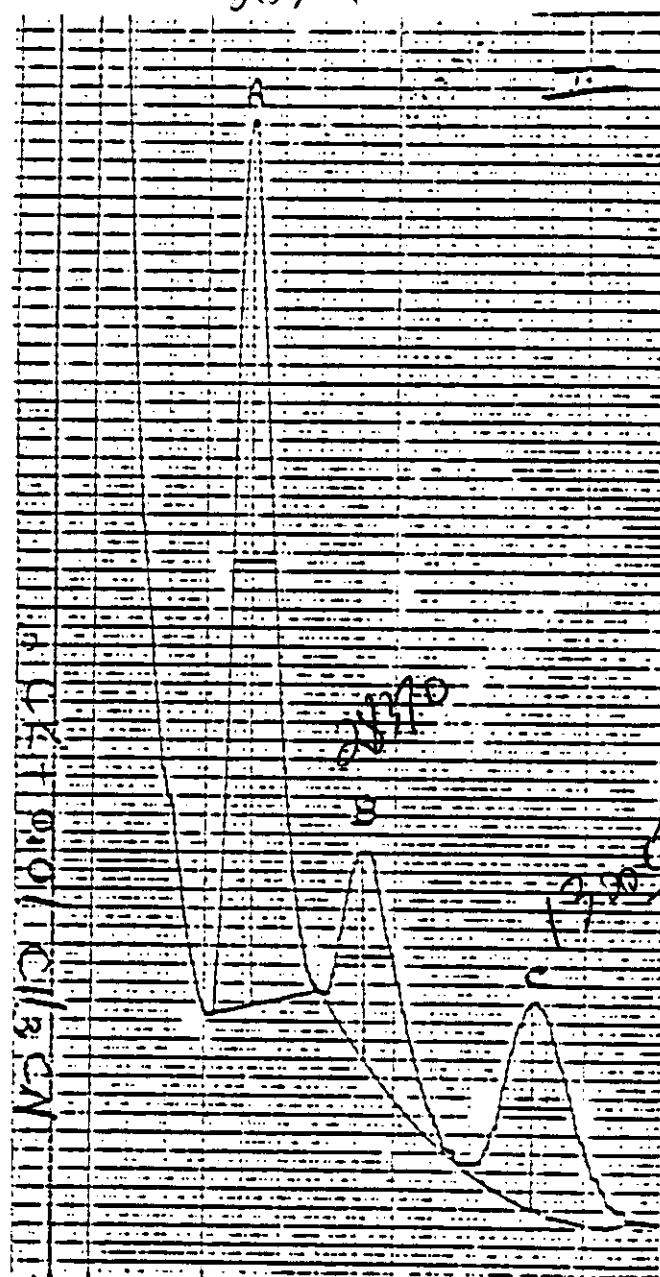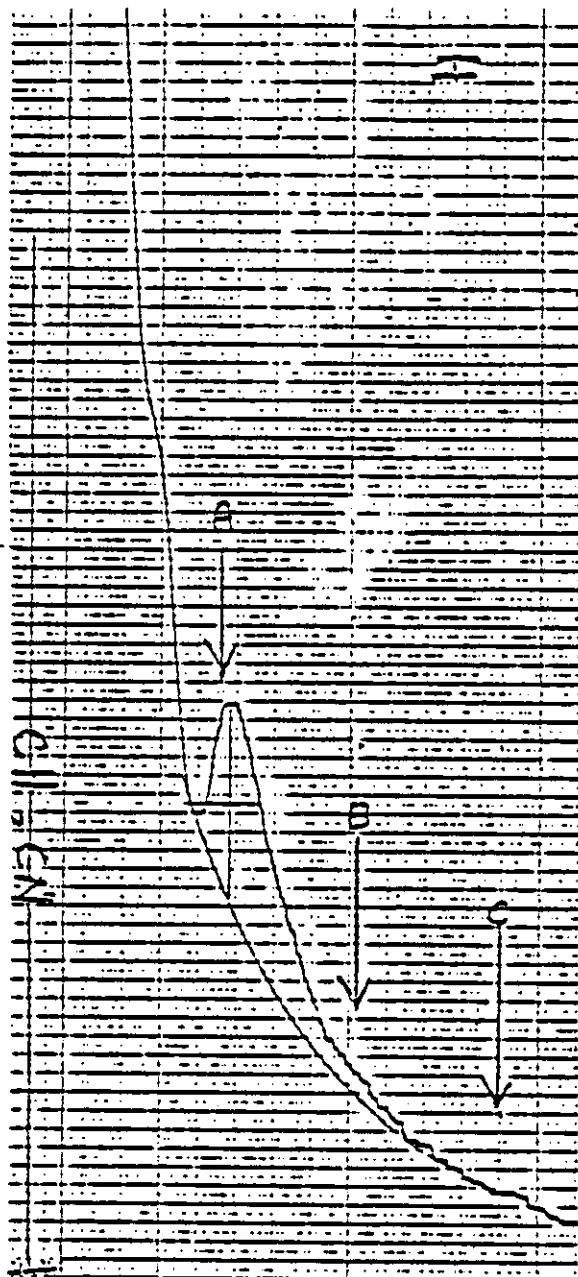

A. GS-28369

B. GS-28370

C. Methidathion

I Check milk CM + 17A

II Check milk plus 0.01 ppm each of the above compounds.

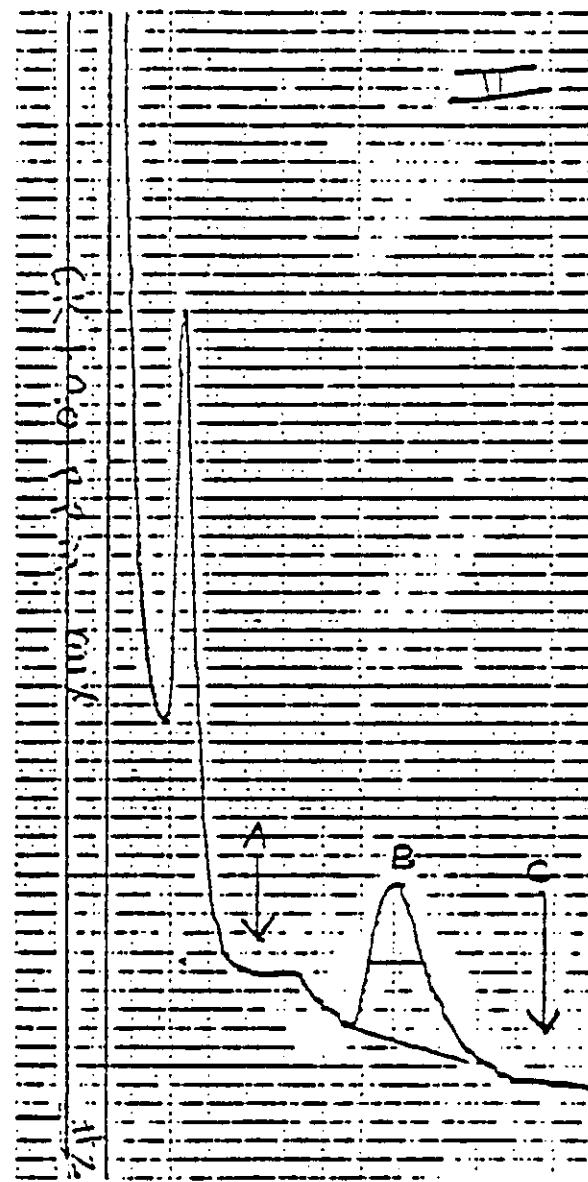
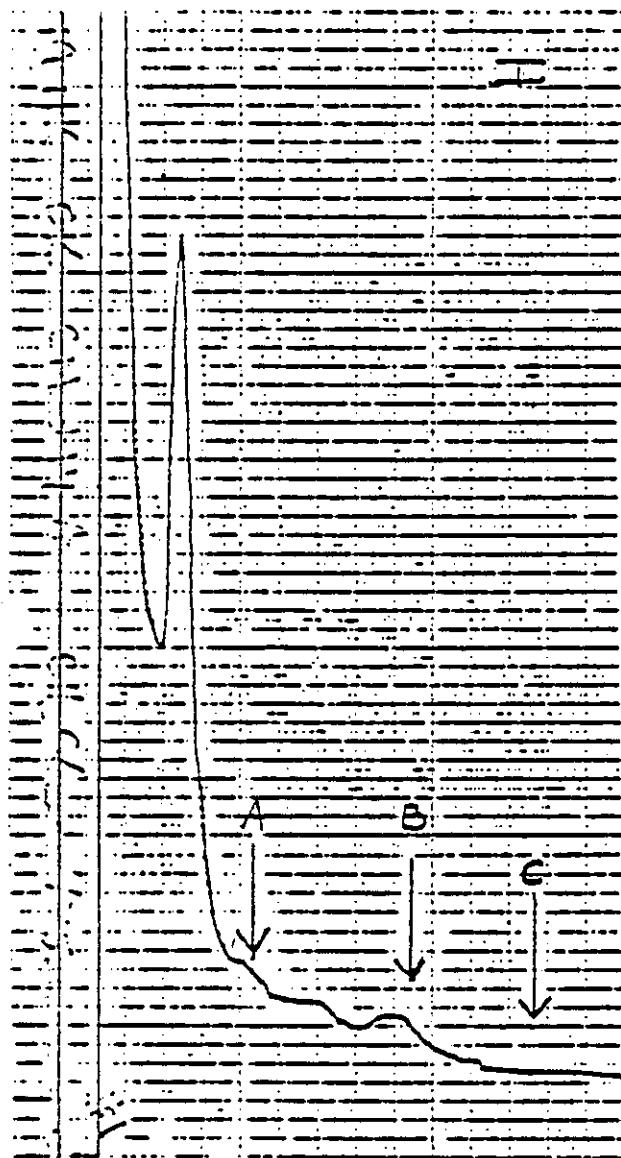


Fig. 6

Gas Chromatograms of Compounds Added to Milk (Methylene Chloride fraction).

- I. Check Milk CM + 17A
- II. Check Milk plus 0.01 ppm of GS-28370.

A. GS-28369 B. GS-28370 C. Methidathion

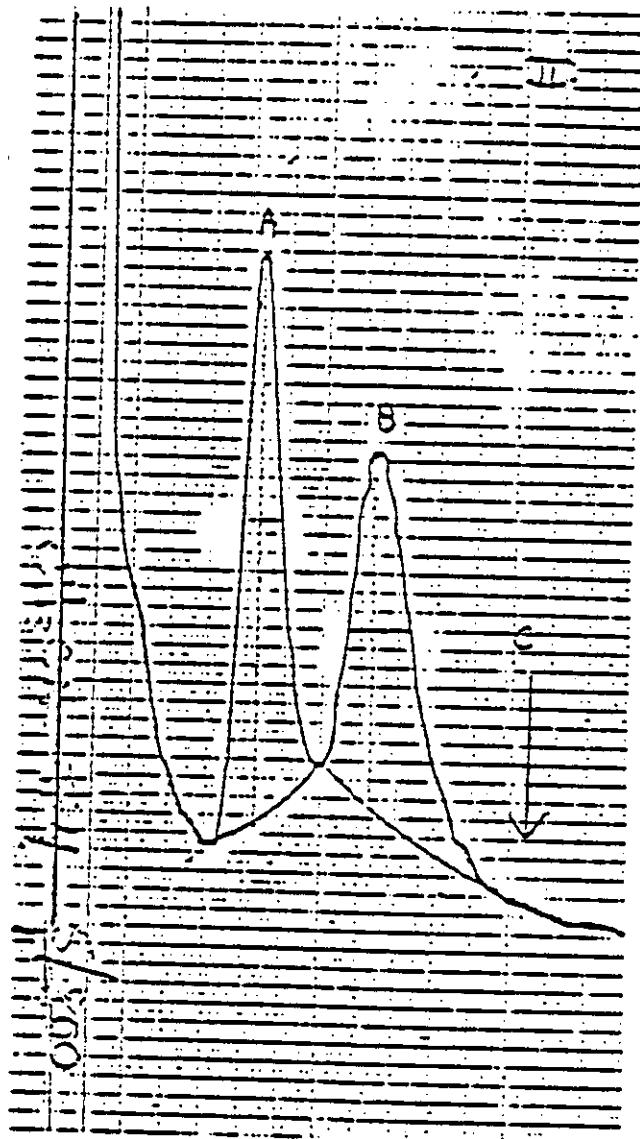
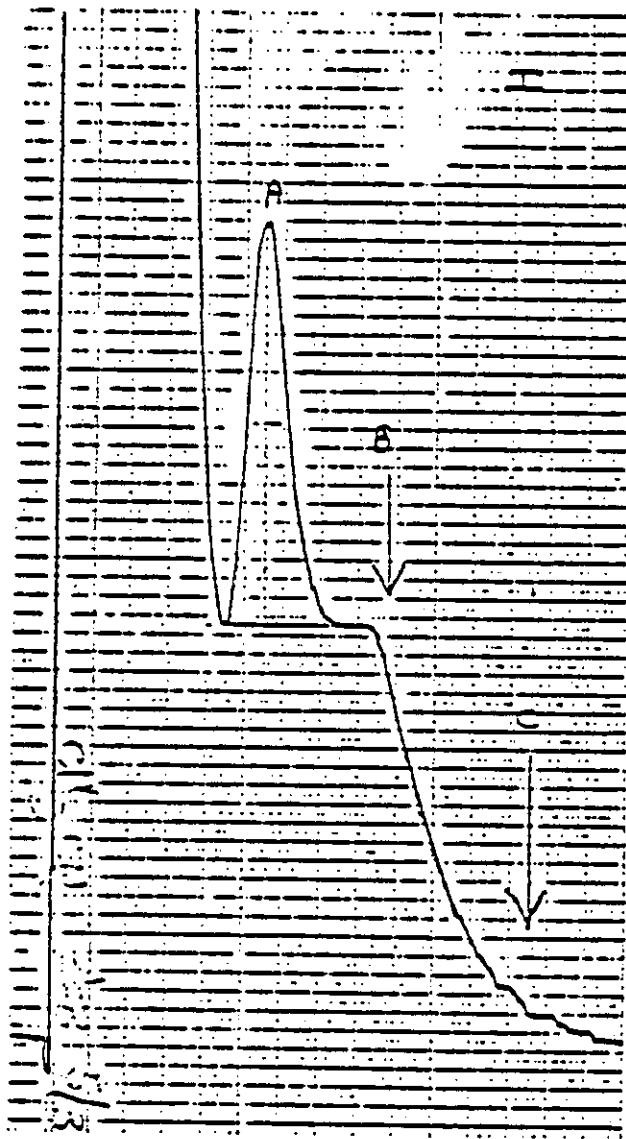


Fig. 7

Gas Chromatograms from Analysis of Treated Milk Samples.  
(Acetonitrile Fraction).

I. Sample Code TM + 17 A  
II. Sample Code TM + 17 P

A. GS-28369

B. GS-28370

C. Methidathion

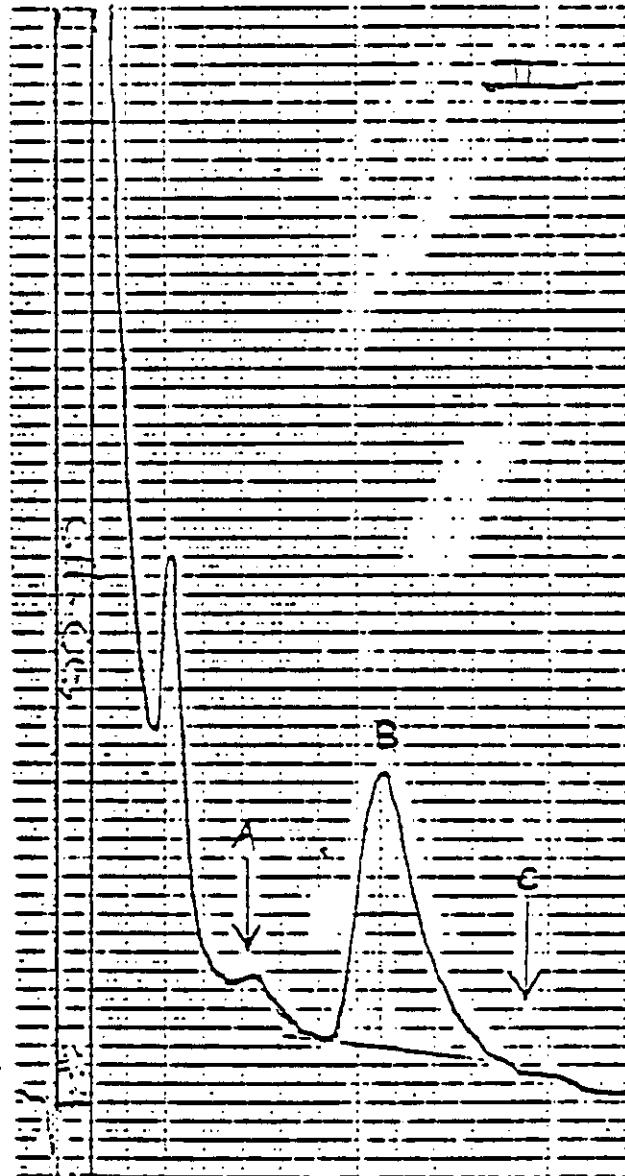
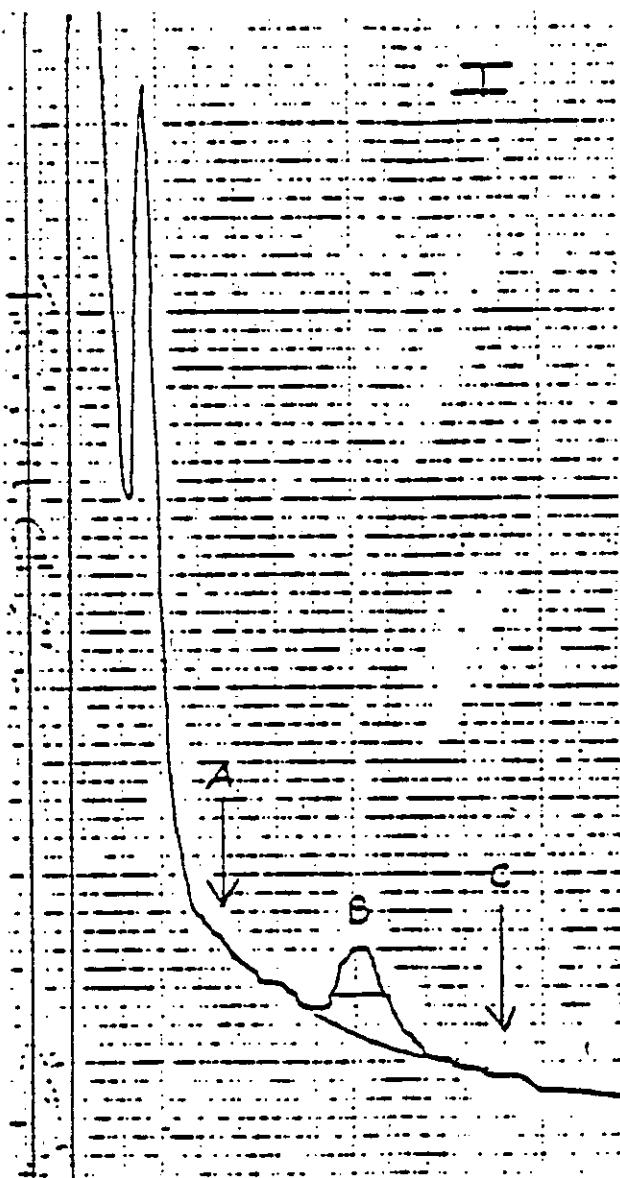


Fig. 8

Gas Chromatograms from Analyses of Treated Milk Samples  
(Methylene chloride fraction).

- I. Sample code TM + 17A  
II. Sample code TM + 17P

A. GS-28369    3. GS-28370    C. Methidathion

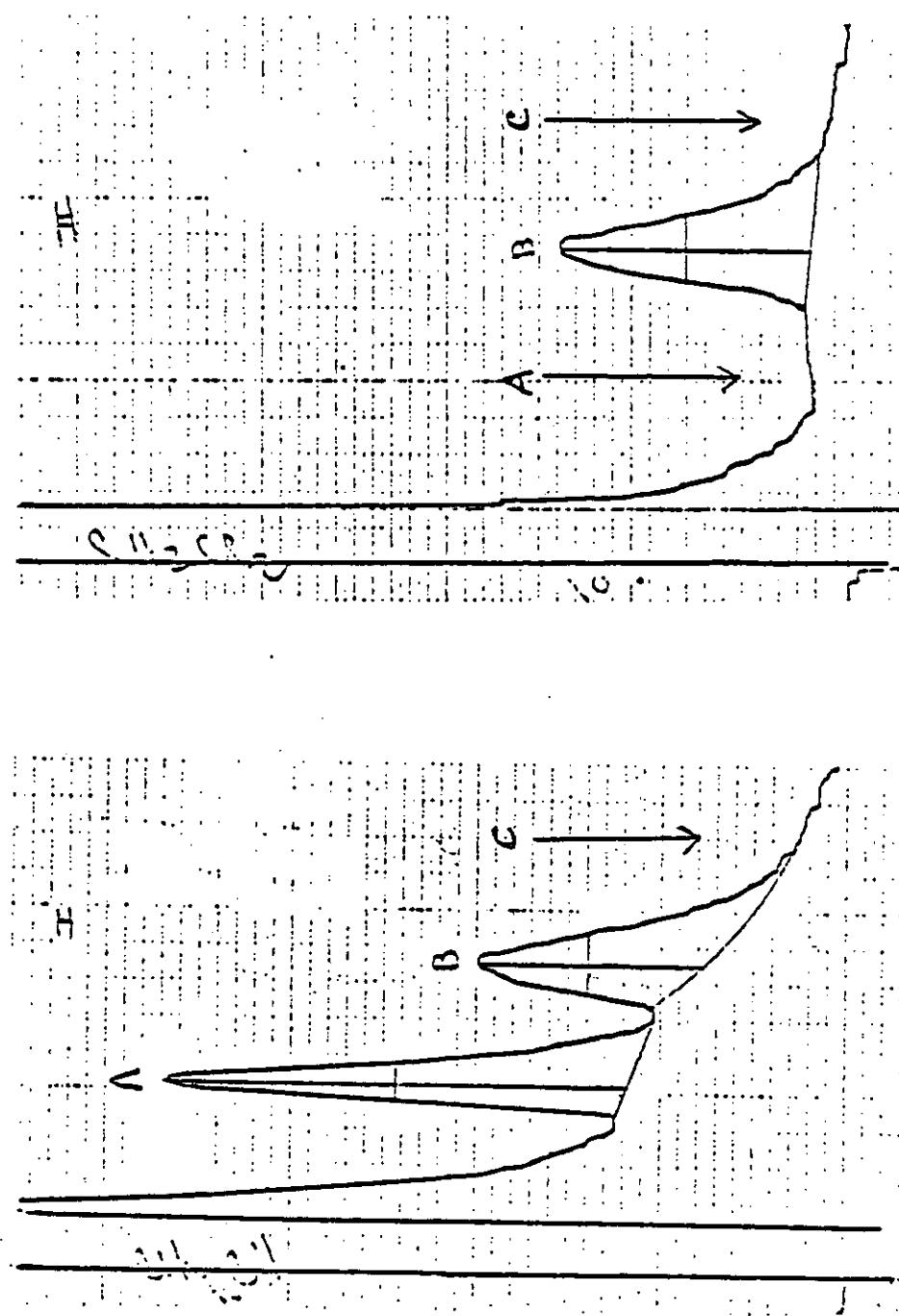


Fig. 9

Gas Chromatograms from Analyses of Treated Milk Samples

Sample Code TM + 25 P

I. Acetonitrile fraction.

II. Methylene chloride fraction

A. GS-28369      B. GS-28370      C. Methidathion

C. Methidathion

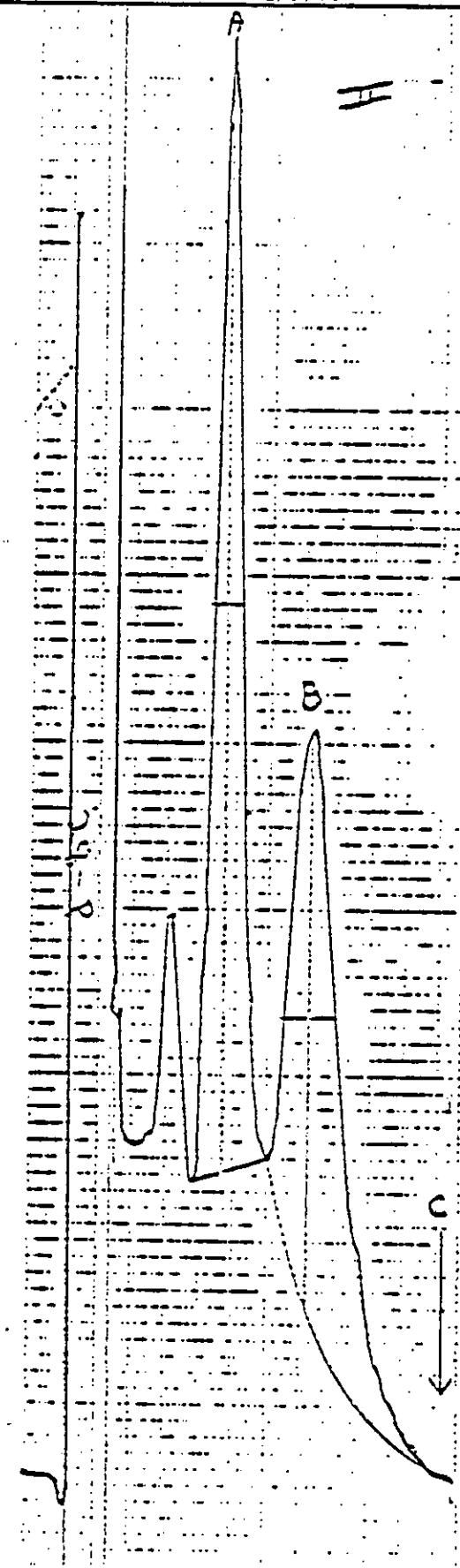
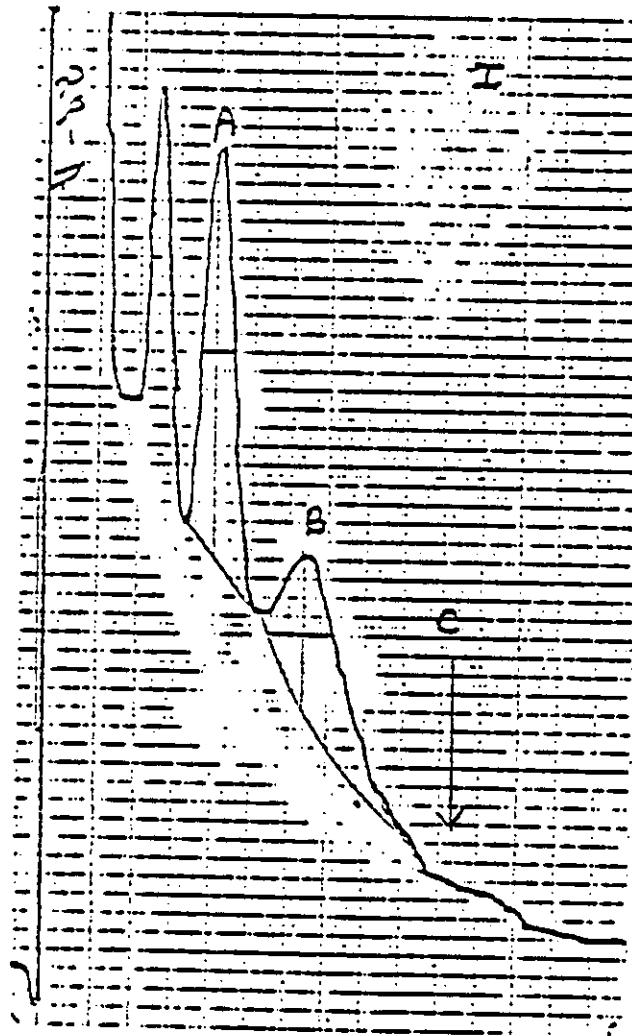


Fig. 10

Gas Chromatograms from Analysis of  
Treated Milk Samples (Acetonitrile  
fraction).

- I. Sample code TM + 23 A
- II. Sample code TM + 23 ?
- A. GS-23869    B. GS-28370
- C. Methidathion

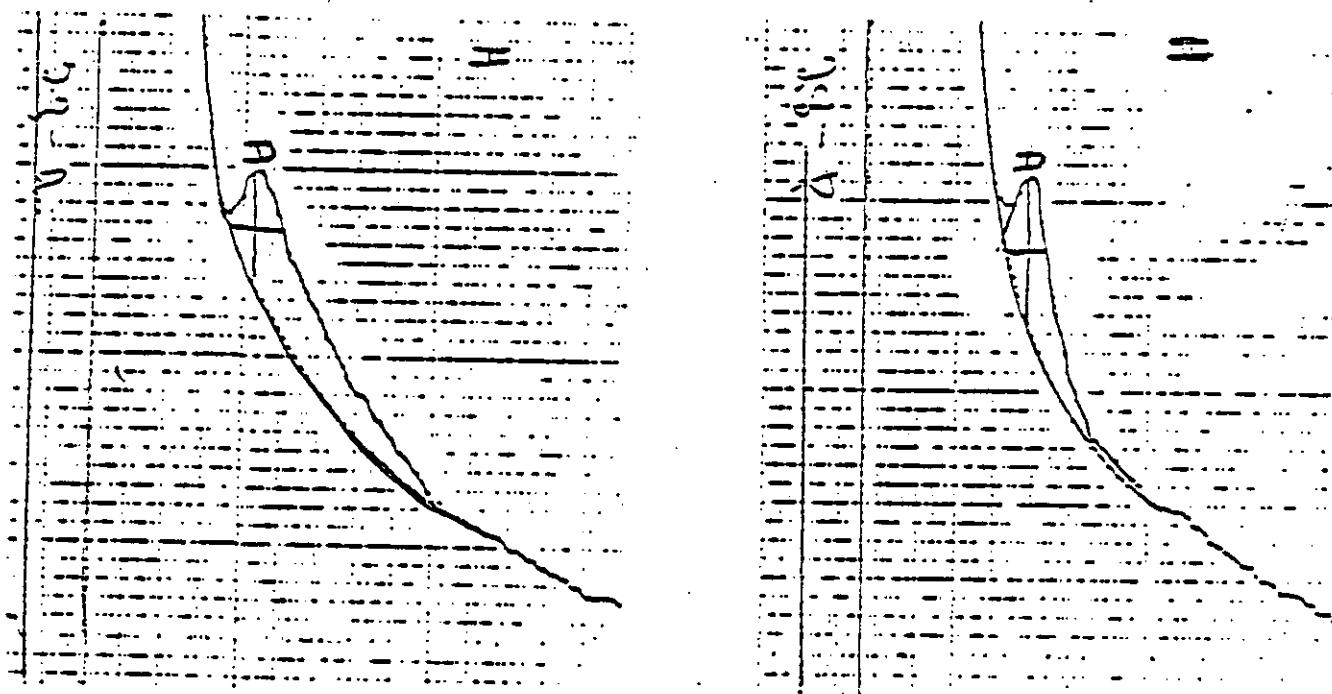


Fig. II

Gas Chromatograms from Analysis of Treated Milk. (Acetonitrile fraction,  
120°C)

- I Sample code TM + 2SA
- II Sample code TM + 2SP

D. GS-28368

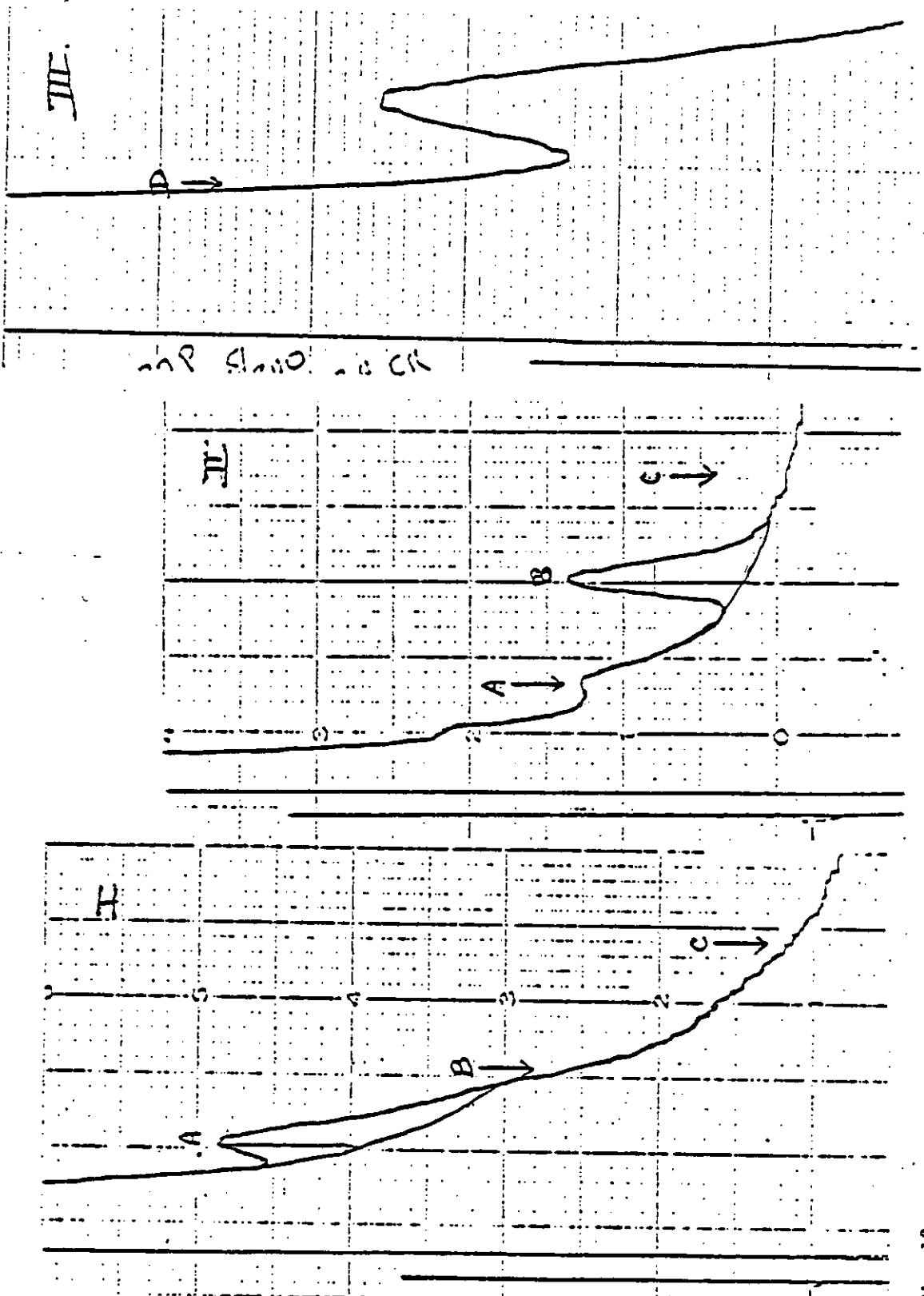


Fig. 12 Gas Chromatograms From Analyses of Treated Milk Samples.


Sample M + 32 P

I. Acetonitrile Fraction

II. Methylene Chloride Fraction

III Acetonitrile Fraction at 120° column temperature.

A. GS-28169      B. GS-28170      C. Methidathion      D. GS-28168

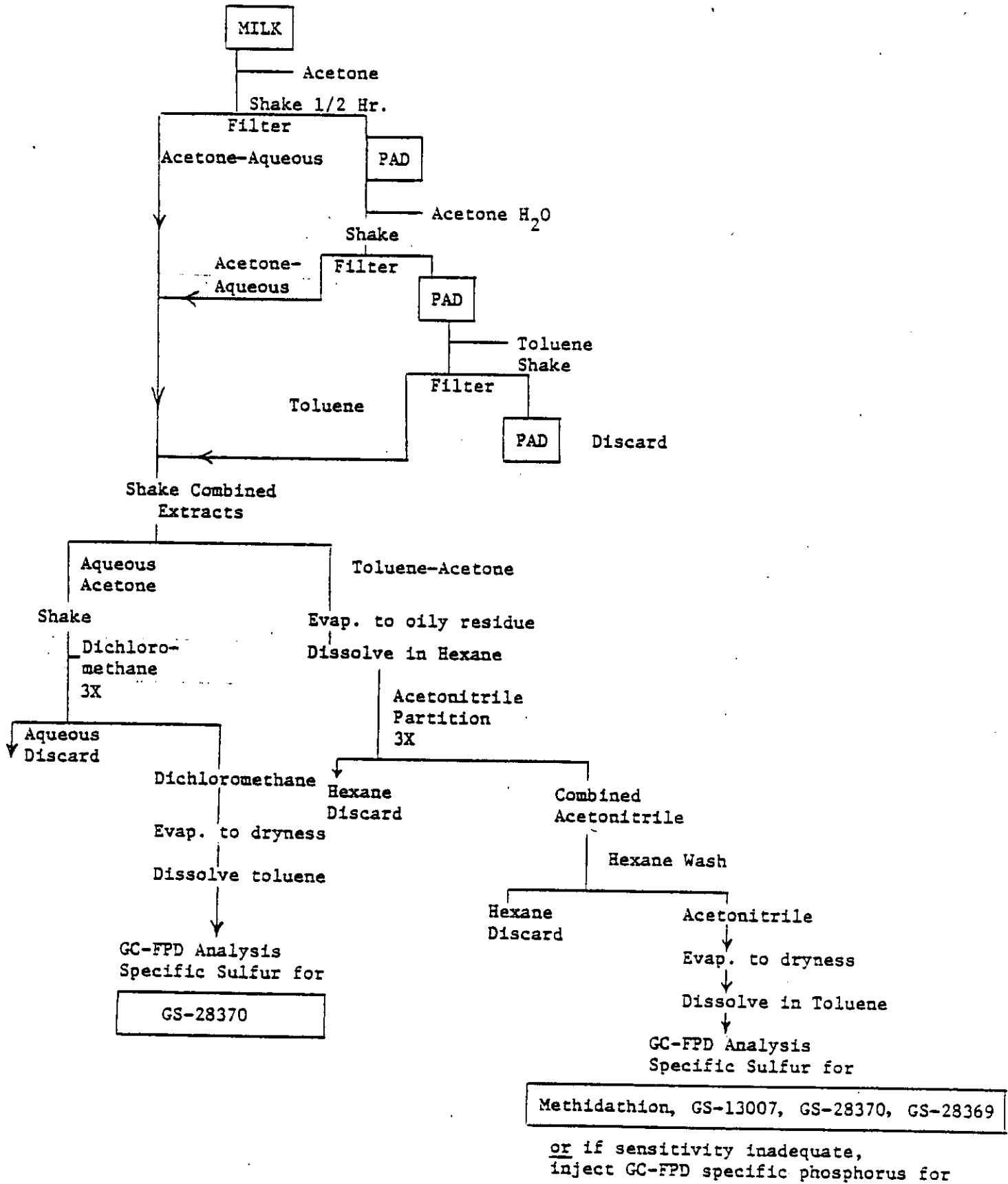

## ADDENDUM TO AG-335

## GC CONDITIONS USED BY ADC LABS

7/15/82

*RAKahn*

All samples were analyzed according to methodology supplied by Ciba-Geigy Corporation. Milk samples were analyzed as described in Method No. AG-335 "Determination of Methidathion and Some of Its Metabolites in Milk by Gas Chromatography". Tissue samples were analyzed according to Method No. AG-334 "Determination of Residues of Methidathion and Its Metabolites, GS-13007, GS-38269, and GS-28370 in Animal Tissue by Gas Chromatography Employing Flame Photometric Detection". These methods were used without modification of any steps, except for filtration of liver samples. Liver extracts were filtered through two pieces of filter paper instead of one.

All gas chromatography was done on a Tracor Model 222 gas chromatograph operated as follows:

Column: 4' x 1/4" i.d. glass packed with 10% DC-200 on Gas Chrom Q, 80/100 mesh

Injection Port Temperature: 239°C  
Column Temperature: 200°C  
Detector Temperature: 219°C

N<sub>2</sub> Carrier Flow (Rotameter): 10  
O<sub>2</sub> Flow (Rotameter): 8  
Air Flow (Rotameter): 100+  
H<sub>2</sub> Flow (Rotameter): 100

Detector Filter: Sulfur for milk; Phosphorous for tissues

Attenuation: 10<sup>3</sup> x 32 for milk; 10<sup>3</sup> x 128 & 10<sup>3</sup> x 64 for tissues

Injection Volume: 5 µl

Chart Speed: 0.25 in/min

Retention Time: GS-13007 - 4.25 min  
GS-13005 - 5.83 min

## ADDENDUM TO AG-335

## GC CONDITIONS USED BY ADC LABS

7/15/82

*RAK/hrs*

All samples were analyzed according to methodology supplied by Ciba-Geigy Corporation. Milk samples were analyzed as described in Method No. AG-335 "Determination of Methidathion and Some of Its Metabolites in Milk by Gas Chromatography". Tissue samples were analyzed according to Method No. AG-334 "Determination of Residues of Methidathion and Its Metabolites, GS-13007, GS-28369, and GS-28370 in Animal Tissue by Gas Chromatography Employing Flame Photometric Detection". These methods were used without modification of any steps.

All gas chromatography was done on a Tracor Model 222 gas chromatograph operated as follows:

Column: 4' x 1/4" i.d. glass packed with 10% DC-200 on  
Gas Chrom Q, 80/100 mesh

Injection Port Temperature: 239°C  
Column Temperature: 163°C  
Detector Temperature: 219°C

N<sub>2</sub> Carrier Flow (Rotameter): 10  
O<sub>2</sub> Flow (Rotameter): 8  
Air Flow (Rotameter): 100+  
H<sub>2</sub> Flow (Rotameter): 100

Detector Filter: Sulfur  
Attenuation: 10<sup>3</sup> x 32 and 10<sup>3</sup> x 64  
Injection Volume: 5 µl  
Chart Speed: 0.25 in/min  
Retention Time: GS-28369 - 5.24 min  
GS-28370 - 6.03 min